


MMA Alarm Cable

Features

- Used to connect medical monitoring devices to the teleCARE IP system via the Multi Medical Alarm module (NIMA)
- Kevlar reinforced cable.
- Fitted with a 2.5mm TS jack plug
- Color gray

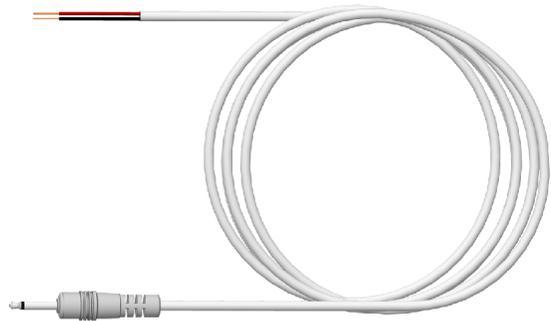

Description

MMA alarm cable with 2.5mm TS jack plug and 4 metres of flexible Kevlar reinforced cable. An article number label is fixed on the cable, below the plug.

Application

This MMA alarm cable is used to connect medical monitoring devices to the teleCARE IP system via the multi medical alarm module (NIMA) and to connect non-medical devices to the multi general alarm module (NIGA).

The two wires (red and black) are used to connect the medical monitoring device. The wires should only be connected to a normally open, potential-free, alarm contact in the medical monitoring device.

In order to fulfil the requirements of the EN 60601-1-8, when used with the NIMA module, an end of line (EOL) resistor must be connected across the alarm contacts of the medical device to monitor the continuity of the electrical connection to the NIMA and the NIGA modules.

Technical Specifications

Markets:		EU	US	UK	ANZ
Item number:	660393	•	•	•	•
Color	Gray				
Length	4 meters				
Conection wires:	red = center				
	black = outer				



Ascom (Sweden) AB
 Grimbodalen 2
 SE-417 49 Göteborg
 Sweden
<http://www.ascom.com/ws>